ULYS5


ULYSS multi-view sinuscope, 10° to 100° 🚈

Ref. 161 401 760

- ➤ Ø 4.2 mm
- ➤ Working length 175 mm
- > 340° shaft rotation
- Autoclavable
- ➤ Titanium body



Plastic sterilization tube for endoscope Ref. 191 402 060



Accessories (in option)



Suction irrigation handle Ref. 191 500 110

Ref. 191 400 110





Plastic sterilization box Ref. 391 308 000


- > For steam sterilization
- > Size 380 x 220 x 80 mm



French production plant: SOPRO-COMEG ZAC Athélia IV Avenue des Genévriers 13705 La Ciotat cedex - France Tél. +33 (0) 442 98 01 01 Fax +33 (0) 442 71 76 90




German production plant: SOPRO-COMEG Dornierstr.55 D-78532 Tuttlingen Germany Tel: +49 (0) 7461-9344-0 Fax: +49 (0) 7461-9344-20

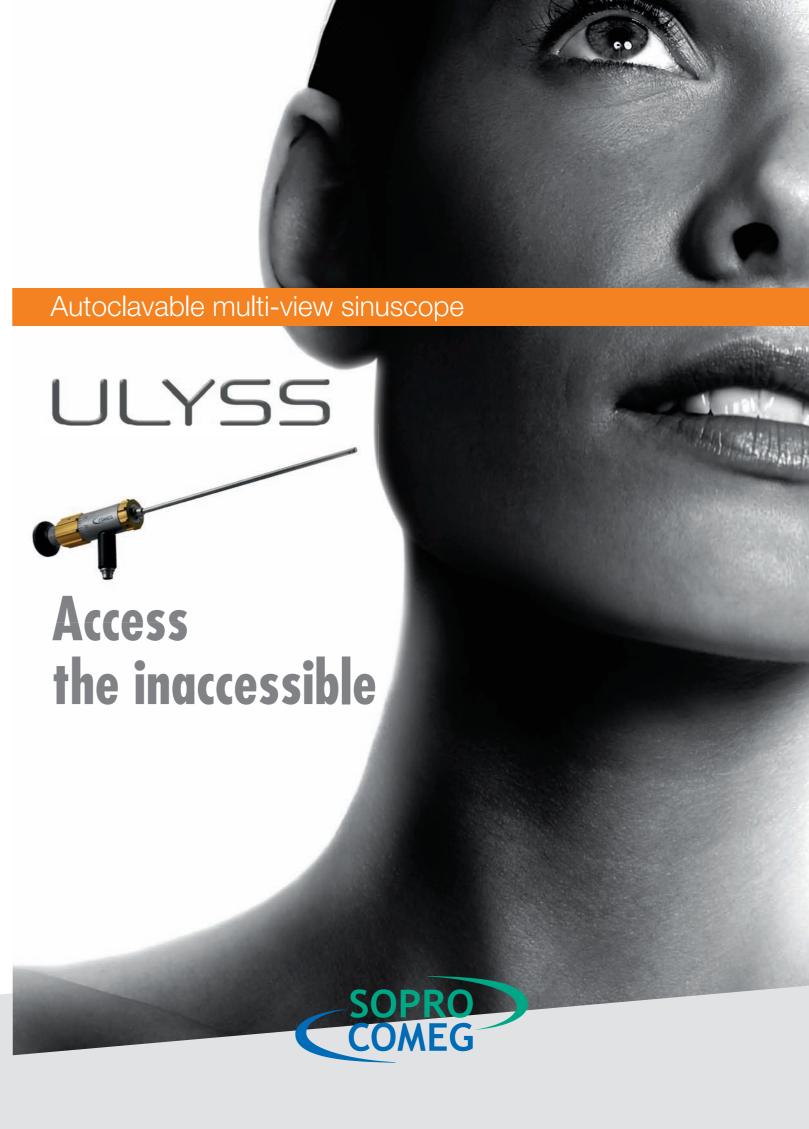



Autoclavable multi-view sinuscope

Access the inaccessible

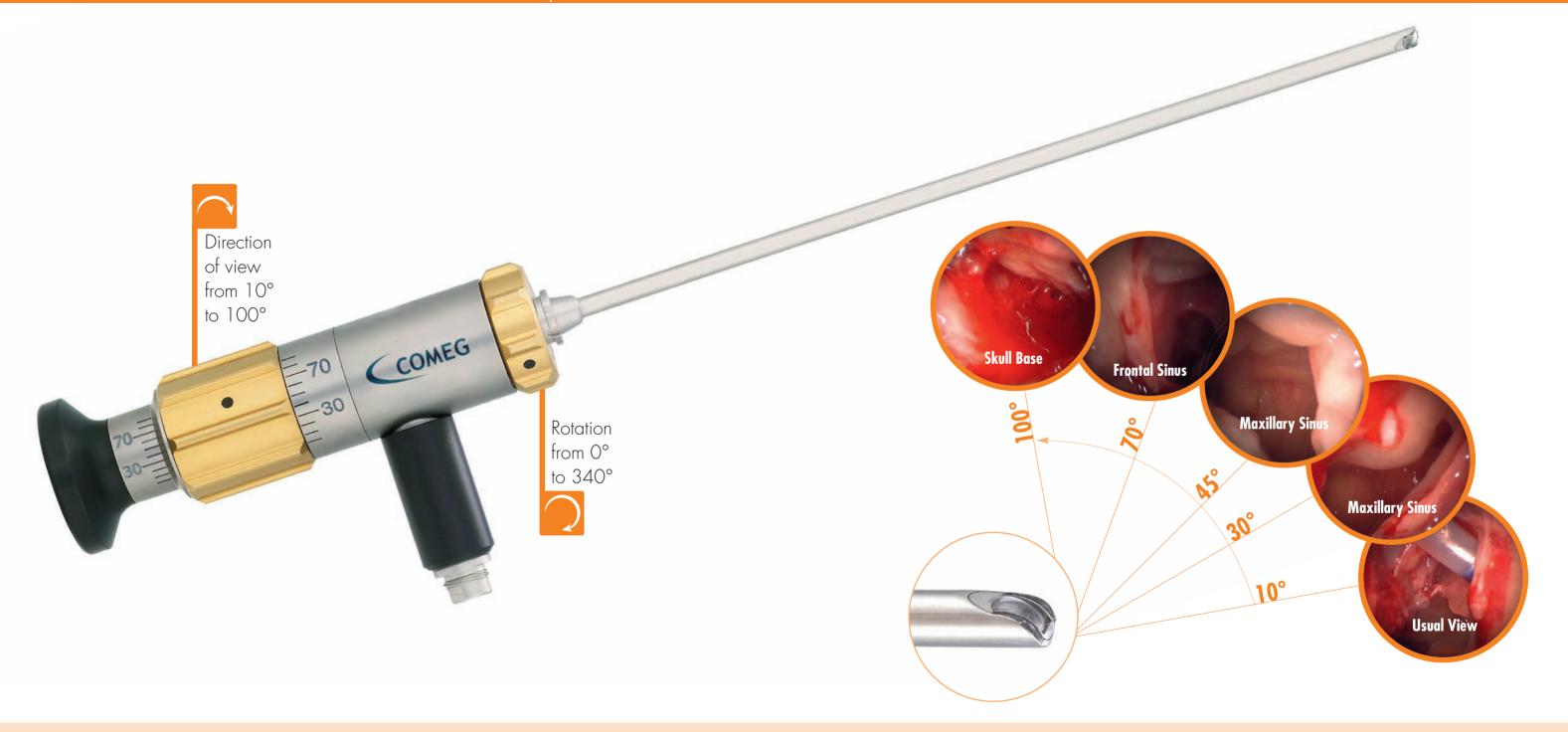

- > ULYSS multi-view sinuscope is the unique autoclavable telescope which allows vizualization of the entire nasal anatomy, including the paranasal sinus spaces. A full diagnosis can be made with the use of only one sinuscope!
- > ULYSS is the perfect solution for surgeries that usually require the use of one or more sinuscopes. The controlled and gradual variation of direction of view facilitates the adaptation of the surgeon to different viewing angles. This allows surgeries such as maxillary sinus mucoceles, balloon sinusplasty, anterior ethmoidectomy and radical sphenoidectomy.
- > The 340° axial rotation associated to the 10° to 100° variable direction of view allows a complete visualization of endo nasal cavity during diagnosis, procedure and verification of surgery with perfect ergonomics.